


ATTO 647N Conjugated GOAT Anti-MOUSE IgG (H&L) Antibody

Product Details

Size	100ul	
Label	iFluor 647	
Species Reactivity	Mouse	
Host/Isotype	Goat / IgG	
Class	Polyclonal	
Туре	Secondary Antibody	
Conjugate	iFluor 647, Ex: 644nm; Em: 669nm.	
Immunogen	Mouse IgG (H+L)	
Form	Liquid	
Purification	Protein affinity purified.	
Storage buffer	Preservative: 0.02% Sodium azide Constituents: 30% Glycerol, 1% BSA, 68.98% PBS	
Contains	0.05% sodium azide	
Storage conditions	Store at +4 $^\circ\!\!\!C$ after thawing. Aliquot store at -20 $^\circ\!\!\!C$. Avoid repeated freeze / thaw cycles	


Applications	Tested Dilution	Publications
Immunocytochemistry (ICC/IF)	1:500-1:1000	-
Flow Cytometry (Flow)	1:500-1:1000	-

Images

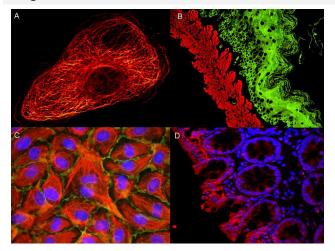

ATTO ® dyes can be used for multicolor immunofluorescent detection with low background and high signal. Examples shown are: A. Tubulin in PtK2- male Rat Kangaroo Kidney Epithelial Cells was detected using ATTO 532 labeled secondary antibody. B. Muscle alpha-actin was stained with a mouse primary antibody and ATTO 488 anti-mouse IgG (green) while Cytokeratin was stained with polyclonal rabbit anti-cytokeratin and ATTO 647N anti-rabbit IgG (red). C. HUVEC (Human umbilical vein endothelial cells were stained with anti- Vimentin-ATTO 532 (green), anti-E-Cadherin-ATTO 655 (red) and DAPI (blue). D. Rat colon sections were stained with Anti-Aquaporin 3-ATTO 594 antibody. Hoechst 33342 (blue) is used as counterstain.

Note: the product information may be optimized and upgraded. Please refer to the actual label information.

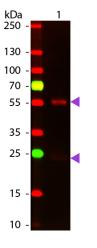

Western Blot of ATTO 647N conjugated Goat anti-Mouse IgG Pre-Adsorbed secondary antibody. Lane 1: Mouse IgG. Lane 2: none. Load: 50 ng per lane. Primary antibody: none. Secondary antibody: ATTO 647N goat secondary antibody at 1:1,000 for 60 min at RT. Block: MB-070 for 30 min at RT. Predicted/Observed size: 25 & 55 kDa, 25 & 55 kDa for Mouse IgG. Other band(s): none.

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will condord or sample node or sample Aurnished to Buyer is and/or accompanying package inserts ("Documentation, the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample Aurnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will condord or sample and/or buyer is metrely ESS OR INPLEOS OR RANNELS OF RENERCHARTISES FOR AMPLEO WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation, accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthori

Note: the product information may be optimized and upgraded. Please refer to the actual label information.